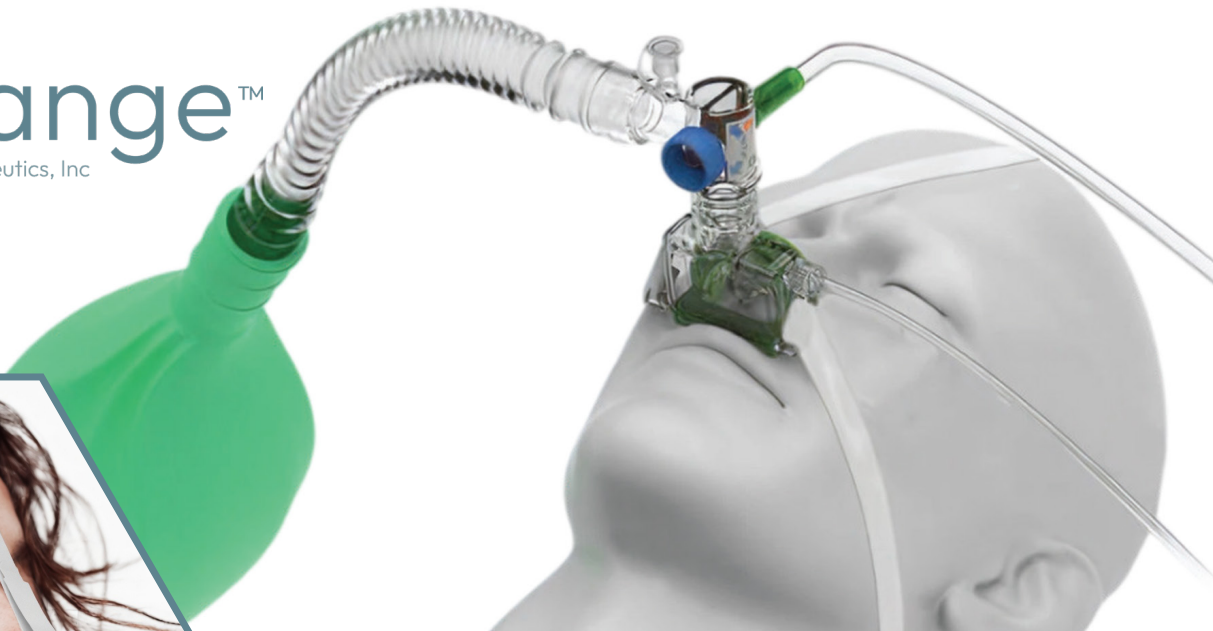




For use in patients with
**High BMI or OSA while under
moderate to deep sedation**

Adjust pressure to pneumatically stent open the upper airway for
superior oxygenation and ventilation

Simple to Use - 3 Easy Steps

- 1 Mount nasal dock and adjust comfort straps
- 2 Connect O₂ line to wall oxygen
- 3 Connect ETCO₂ line to monitor

*Now you can adjust and confirm pressures with the integrated APL
valve and manometer

Feature	xchange™	Standard Nasal Cannula	High-Flow NC
Max CPAP Pressure	20+ cmH ₂ O	Limited	Variable 1-7cmH ₂ O
Visually Confirm Pressure	Yes - manometer	No	No
EtCO ₂ Sampling	Consistent	Often unreliable	Variable
Requires Capital Expense	No	No	Yes
Patient Comfort (Nasal)	High	Moderate	Moderate

